

Cap2v8® Sterile Media Bottle Assembly



Engineered, Single-Use Solutions —
Reducing Operational Cost, Improving
Production Efficiency

- Molded and Unitized Construction (No Adhesives)
- Certified Sterile with Fast Delivery
- Cost-Effective Aseptic Transfer and Sampling Device
- One Continuous Flow Path
- Lot Traceable
- Meets USP Class VI

FEATURES

BENEFITS

One-Piece Molded Cap Design

Reduces contamination and product loss

Materials of Construction

C-Flex® tubing (heat weldable/sealable); Polypropylene cap; PETG Bottle


Single Use

Cost effective; Reduced assembly and installation time; Eliminates cleaning validation

Custom Design Capabilities

SaniSure® offers a full array of custom options, including tubing, fittings/end connections, sterile connectors and bottle variations

SaniSure offers a full array of custom options, including tubing, fittings/end connections, sterile connectors and bottle variations.

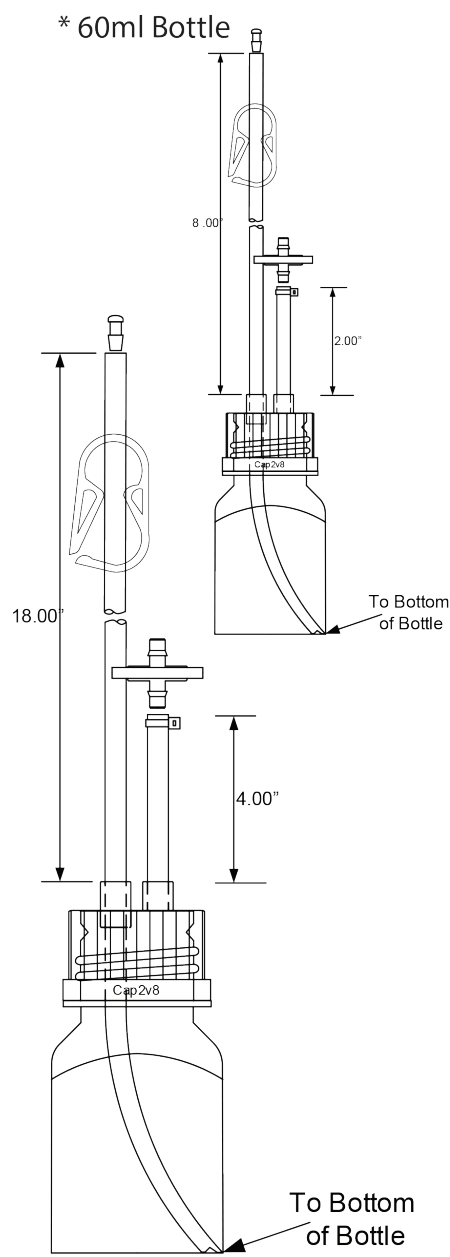

Packaging

Bagged and Labeled in an ISO Class 7 Cleanroom to Your Specifications.

Applications


- Tissue and cell culture media
- Cell fermentation
- Pilot and process research and development
- Laboratory research
- Aseptic sampling
- Vaccine production
- Sterile filling and transfer

Sterilization

Gamma irradiated to sterility assurance level 10^{-6}

| Part Number | Bottle Volume | Bottle Material |
|---------------------|---------------|-----------------|
| CAP2V8-60PETG-X-2 | 60 mL | Square PETG |
| CAP2V8-125PETG-X-2 | 125 mL | Square PETG |
| CAP2V8-250PETG-X-2 | 250 mL | Square PETG |
| CAP2V8-500PETG-X-2 | 500 mL | Square PETG |
| CAP2V8-1000PETG-X-2 | 1 Liter | Square PETG |
| CAP2V8-2000PC-X-2 | 2 Liter | Round PC |

Bottle assembly includes: 25 mm 0.2µm Vent Filter, Pinch Clamp and Tubing Plug, transfer line tubing V-notched to reduce/eliminate suction to bottom of bottle

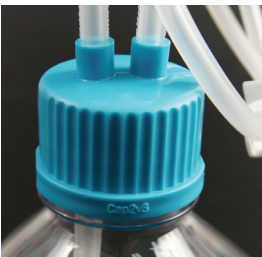

SaniSure Cap2v8® is now available with AdvantaFlex® — custom colored cap for easy identification.



IMPORTANT: It is the user's responsibility to ensure the suitability and safety of SaniSure, Inc. products for all intended uses.

For a period of 90 days from the date of delivery, SaniSure, Inc. warrants this product to be free from defects in materials and workmanship. Our only obligation will be to replace any portion proving defective, or at our option, refund the purchase price thereof. User assumes all other risk, if any, including the risk of injury, loss or damage, direct or consequential, arising out of the use, misuse or inability to use this product. THIS LIMITED WARRANTY IS IN LIEU OF THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR PARTICULAR PURPOSE AND ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED. No deviation is authorized.

SaniSure, Inc. assumes no obligations or liability for any advice furnished by it or for results obtained with respect to those products. All such advice is given and accepted at the user's risk. SaniSure® is a regis-